



User's Quick Install Guide

Instructions for Installing Your New Personal Emergency Response System

Instructions:

1. Unpack the LifeSentry from the box and confirm that you have all of the pieces and accessories. (See Step 1 of this guide.)
2. Follow the step-by-step instructions on the following page to set up your LifeSentry. (See Steps 2-3 of the guide.)
3. **Put on your pendant.** If you do not wear it, you cannot press the button to reach help in an emergency.
4. Press the blue emergency button on the front of your pendant to test the alarm system. (See Step 4 of this guide.)

NOTE: When you are testing your LifeSentry, tell the operator your name and that you are a new client.

NOTE: If the test is not successful, contact your dealer.

ATTENTION

- Be sure to fully charge one of the white pendant batteries prior to placing your first test call to the central station. (See step 2 of this guide.)
- Place the fully charged battery in the pendent. (See step 2 of this guide.)
- We recommend that you test the battery and system regularly by pressing the gray battery test button on the back of the pendant. (See Step 4 of this guide.)
- Place the spare pendant battery in the charger and keep it there for when you need to change batteries. Don't worry, it is ok to leave the battery in the charger until you need it. The charger is designed to keep the battery from overcharging. The battery in the pendant will stand by for many months. We recommend that you exchange the batteries monthly.
- Contact your dealer for warranty information.

Equipment Reference:





User's Step-by-Step Installation



1

Confirm Accessories in Kit

4 Rechargeable AA Backup Batteries



2 Lithium Ion Pendant batteries



Telephone Cord



AC Adapter



Lanyard



Wrist strap



Beltclip



2

Place Batteries in Unit



Place Green Rechargeable batteries in the compartment on the bottom of the base unit. **Make sure to turn ON/OFF switch to ON. You will hear the announcement "Running on Battery Power." Unit will do so until the AC Adapter power supply is plugged in.**



Place spare pendant Lithium Ion battery in the charger on the front of the base unit.

(Make sure charging light turns RED, indicating that the charger is ON and that the battery is being charged.)



Place one Lithium Ion battery into the pendant by unscrewing the cap and inserting FLAT SIDE UP.


As soon as cap is tightened, the red light flashes as the pendant does a self-check and announces the battery status.

3

Plug in the Telephone and Power Cords



Plug the AC Adapter power cord into the power jack and twist right to lock.

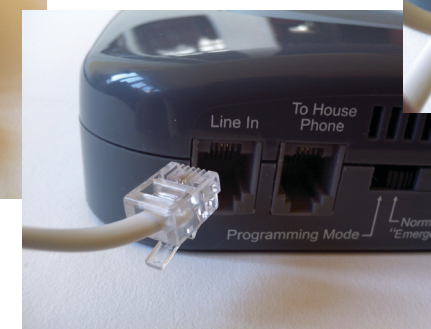

Plug one end of the supplied phone cord into the "Line In" jack. Connect the other end into the wall (or modem) phone jack.



(OPTIONAL) Plug the telephone into the the "To House Phone" jack.

4

Test the Unit

Place a call to the Central Station

1. Press the Blue Button for 4 seconds to send test alarm
2. The Red light will flash. You will hear dialing and other sounds as the alarm is sent.
3. When the operator answers, just tell them you are testing.
4. Once you speak to the operator, setup is complete!



Perform Bi-weekly system tests

Press the gray Test Button for 4 seconds, or until you hear two beeps. This test informs of base/pendant communication, telephone line detection, and battery status.

Perform this test every other week or if the INFO light is blinking.

